

PLUG & TOUCH BY RODENSTOCK

CHART PROJECTOR
RODACHART 360

Displays efficiency to the max.



■ Energy-efficient LED illumination

■ Brighter and clearer chart display

■ Sophisticated design

■ Quick chart selection

■ Programmable display sequence

QUALITY IN DETAIL








Modern usability design

The Rodenstock LED **Chart Projector Rodachart 360** is designed for durability and efficient operation.

An energy-efficient LED light source with 500,000 hours lamp life ensures maintenance free operation with the added advantage of a brighter, high-contrast display with minimum running costs.

The intuitive remote control helpfully allows the programming of three chart sequences – each with up to individual charts – chosen from the 36 different charts available.

The **Rodachart 360** comes complete with remote control, projection screen and sturdy base allowing it to be easily set up for immediate operation. Alternatively the projector can be fitted to a refraction unit or wall-mounted.

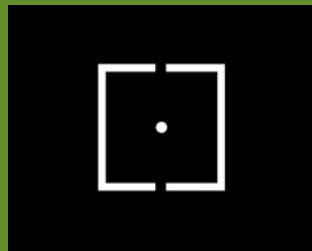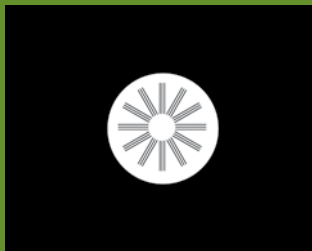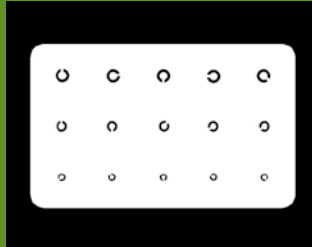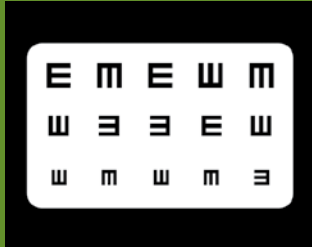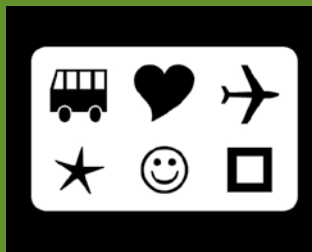

SPECIFICATIONS

Chart change time 0.014 s
 Auto off time 5/10/15 min
 3 user programs 15 steps each
 Masks Single, horizontal, vertical, red/green, polarising
 Projection range..... 2.9 to 6.1 m

Projection

Magnification..... 30x (at 5 m)
 Lamp..... 3 W LED (373 to 801 nm)
 Background luminance 200 cd/m² at 5 m (adjustable)
 Projection tilt Approx ± 10°

Dimensions WDH..... 191 x 310 x 200 cm

Weight..... 4 kg

Power consumption..... 11 W

Input AC 115 V to 230 V

Frequency..... 50/60 Hz

Remote control Included
 (3 channels adjustable)



Remote control

201106 - subject to change without notice