

TR Classic

Tilt table



Specifications

Cushion length - Standard	1750 mm (69")
Cushion length - Pediatric	1400 mm (57")
Cushion width - Standard	650 mm (24 1/2")
Cushion width - Pediatric	450 mm (17 3/4")
Height (Horizontal position)	600 mm (24")
Working Load Limit	150 kg (330 lbs)
Power supply	230 V 50 Hz, 115 V 60 Hz
Power failure lowering	Battery backup (9 V)
Frame	Powder coated steel
Castors	125 mm (5") Two castors with brakes and two with straight steering
Control	Hand control
Cushion	Foam filled, vinyl cover



Registered at Medical Products Agency, Sweden, as a MD Class 1 product, MDD 93/42EEC.

The TR Classic Tilt Table provides a smooth and effortless lift from horizontal to upright position

The TR Classic Tilt Table is designed to make testing, rehabilitation, positioning and transferring patients easier and safer. The TR Classic Tilt Table is a logical choice for patients who have paralysis, are bedridden or wheelchair bound who would benefit from passive standing to improve circulation or reduction of lower extremity edema.

The TR Classic Tilt Table can be used in all environments to assist the caregiver with transfer and treatments needs. The TR Classic Tilt Table operates on a standard 230 V or 115 V power source for smooth and silent operation.

Standard equipment includes three safety belts for comfort and security, footboard, removable wire basket and hand operated control.



TR Classic has a wide foot plate, adjustable and lockable in three angle positions.



TR Classic Tilt Table is operated by an electrical actuator with battery back up for lowering in case of power failure.

TR EQUIPMENT

TR EQUIPMENT AB PO Box 116, S-573 22, TRANÅS, SWEDEN
 TR GROUP INC. 903 WEDEL LANE, GLENVIEW, ILLINOIS 60025, USA
 TR EQUIPMENT UK LTD. INNOVATION CENTRE, GALLOWES HILL, WARWICK CV34 6UW ENGLAND
 TRANSPORHAMED GMBH GARTENSTRASSE 9. 42799 LEICHLINGEN. GERMANY

Phone: +46 140-38 50 60 Fax: +46 140-38 50 61
 Phone: +1 800-752-6900 Fax: +1 847-724-8566
 Phone: +44 844 335 8386 Fax: +44 844 335 8383
 Phone: +49 2175-16 94 22-0 Fax: +49 2175-16 94 22-19

Hygiene Equipment Designed for People

www.trequipment.com



TR Classic

Tilt table



Features

The TR Classic Tilt Table is electrically operated for smooth and silent operation. It effortlessly tilts from horizontal to a 90 degree upright position.

A hand control is provided for convenient tilt adjustment and tilt indicator for precise monitoring of positioning.

The wide footplate, adjustable and lockable in three angle positions, meets the varied needs of the patient population.

The table comes standard with a back-up battery for lowering in case of a power failure. Wire storage basket is removable.

The TR Classic Tilt Table easily manoeuvres on four 125 mm (5") locking castors with straight steering. Wheel fenders protect the lying surface when transporting through narrow doorways.


Models

Standard

Pediatric

Accessories

Activity Tray

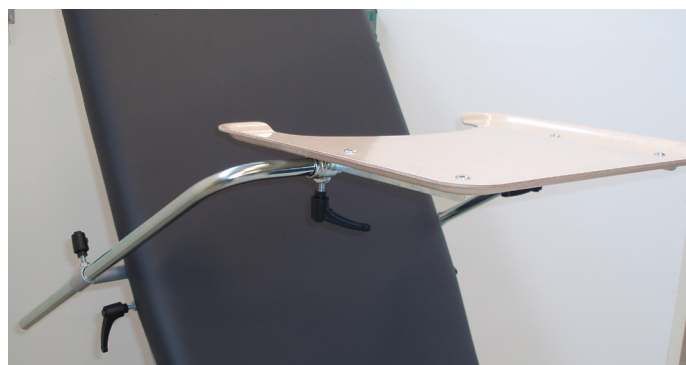

TR Classic effortlessly tilts from horizontal to a 90 degree upright position and has a tilt indicator for monitoring of angle.

TR reserves the right to make design changes and/or improvements without notice

Distributor:

TR EQUIPMENT

TR EQUIPMENT AB

PO Box 116, S-573 22, TRANÅS, SWEDEN

Phone: +46 140-38 50 60

Fax: +46 140-38 50 61

TR GROUP INC.

903 WEDEL LANE, GLENVIEW, ILLINOIS 60025, USA

Phone: +1 800-752-6900

Fax: +1 847-724-8566

TR EQUIPMENT UK LTD.

INNOVATION CENTRE, GALLOW'S HILL, WARWICK CV34 6UW ENGLAND

Phone: +44 844 335 8386

Fax: +44 844 335 8383

TRANSPORHEMED GMBH

GARTENSTRASSE 9, 42799 LEICHLINGEN, GERMANY

Phone: +49 2175-16 94 22-0

Fax: +49 2175-16 94 22-19

Hygiene Equipment Designed for People

www.trequipment.com

