**DC-60 Exp** Diagnostic Ultrasound System

# An all-around partner to trust

X-Insight





www.mindray.com

P/N:ENG-DC-60 Exp-285210X12P-20180706 ©2018 Shenzhen Mindray Bio-Medical Electronics Co.,Ltd. All rights reserved.







## An insightful solution to envision more.


Mindray's brand-new solution is an excellent fusion of continuous customer insights into clinical needs, combined with continually evolving cutting-edge ultrasound technologies.

Full of vitality, with the goal of insight into the future and pursuit of the infinity, it is constantly improving with increased scalability.

## High efficiency with precision imaging

As an all-around partner, DC-60 Exp with X-Insight devotes to provide a comprehensive solution to help you manage all-aspect daily clinical practice with ease and certainty.

Based on the deep insights of customer needs, the DC-60 Exp with X-Insight is designed to deliver high efficiency with precision imaging, which is empowered by eXpress Clarity, eXceptional Intelligence and eXceeding Experience :

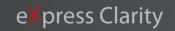



eXceptional Intelligence

eXceeding Experience



# eXpress Clarity


More clarity at hand

To achieve excellent images with minimum effort is your highest priority and our endless pursuit. With continuous innovation in imaging and transducer technologies, the DC-60 Exp with X-Insight delivers immediate clarity with 'Touch and Bingo' to allow you to get the optimum images that you need immediately, as soon as the transducer touches the body.


### 'Touch and Bingo' with evolutional transducers

### ComboWave transducers

Compared with traditional transducers, ComboWave transducers utilize a new type of composite piezoelectric material to dramatically optimize the acoustic spectrum and reduce acoustic impedance. Further integrated with Mindray unique 3T technology, the ComboWave linear transducers allow you to experience outstanding performance with extreme image resolution and uniformity in thyroid, breast, vascular, and more.

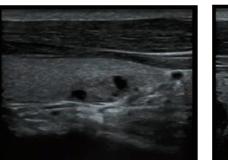

Traditional Linear Image

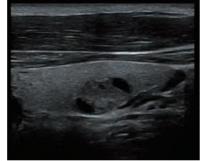

ComboWave Linear Image

## Thunder-speed imaging powered by X-Engine

The new X-Engine utilize the evolutionary integration of GPU and CPU to enables multi-core parallel processing for fast imaging. With the advanced imaging engine, the imaging processing speed is accelerated three or four times faster than traditional processing, resulting in extremely fast imaging for 3D/4D and other applications.

## Ultra-Speed data transmission with security guarantee


By the integration of Solid-State Disk(SSD), DC-60 Exp with X-Insight provides the ultra-speed data transmission capacity to accelerate the system operational speed, for instance, faster imaging, quicker boot-up or wake-up. And the SSD is able to extends the data life-time of ultrasound system, which enhances the protection of patient data with safety and privacy guarantee.




_	
-	
_	
_	

### Optimized transducers with 3T technology


Combining with Mindray unique 3T technology (Triple-matching layers, Total-cut design, Thermal control), the optimized convex, phased array and volume transducers provide a wider full range application coverage, resulting in an optimum scanning solution in ABD, Cardiology, OB/GYN and more.

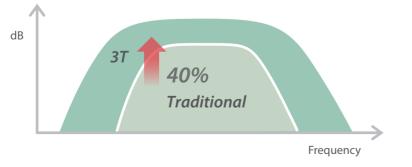

### Proven advanced application:

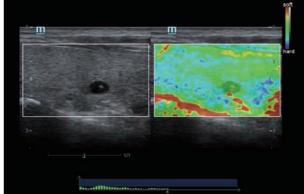

Natural Touch Elastography

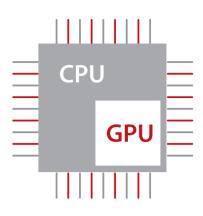

Processing Speed 3x Plus



data

Transmission speed 5x Plus



UWN+(CEUS)

# Image gallery




Liver, 2D



Abdomen Mass



Thyroid, 2D

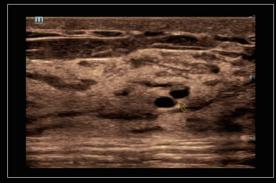


Dilation of Mammary Duct



Kindy HR Flow




Bladder Mass, 2D

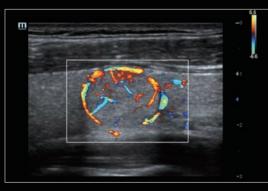


Thyroid Mass CFM



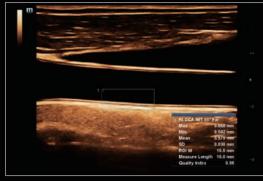

Common Carotid Artery IMT

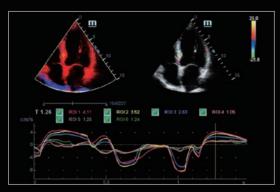

Cardiac TDI QA




Fetal Leg 3D




Fetal Heart iLive

Subscapular Muscle Tendon

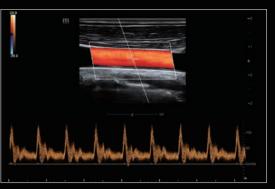


Common Carotid Artery PW

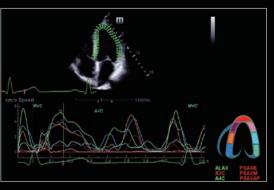

Cardiac TT QA



Fetal Face iLive

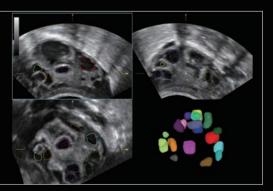

Smart Follicle



## eXceptional Intelligence


### Intelligence through whole exams

To improve scanning efficiency with more accuracy and consistency, the DC-60 Exp with X-Insight provides All-Smart Exams for full kinds of advanced application, from plane acquisition to image optimization, and from calculation to exam protocol.

### **Smart Track**

More simplicity with auto image optimization

Smart Track provides a fast and intelligent optimization for vascular imaging with simply **one-touch operation**. It can optimize Color, Power and PW spectrum by **self-tracking** and reduce time-consuming steps. Therefore the vascular exam workflow is simplified with an optimal view.

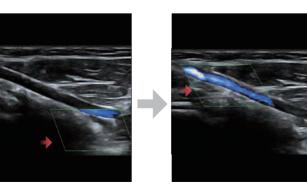

Smart Track

### **Smart Planes CNS**

More accuracy with smart acquisition

Smart Planes CNS is a user-friendly to improve through-put, and reduced user dependency. **With a simple button-click** on a 3D fetal brain volume image, the standard CNS scanning planes (MSP, TCP, TTP and TVP) and a range of related anatomical measurements (BPD, HC, OFD, TCD, CM and LVW) are obtained immediately with high accuracy, which reduces about 90 Second per exam.

### Auto EF More productivity with auto calculation

Auto EF is an intelligent way to analyze 2D echo clips to **automatically recognize** diastole & systole frames and **output** a series of measurements to evaluate left ventricle function for more productivity, such as EDV/ESV/EF.


### iWorks

#### More consistency with exam protocol


iWorks takes advantage of the inbuilt standard scan protocols for more consistency and **reduces the exam time by up to 50%**. The highly flexible and powerful feature is further enhanced with user-defined capabilities.



Smart Planes CNS

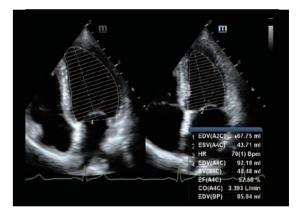

Auto EF

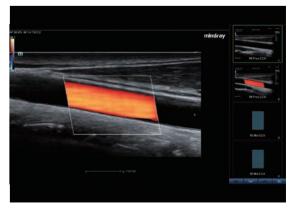

iWorks

## eXceeding Experience

Experience with high productivity

Comfortable experience in scanning results in more patient focus for high quality exams in both clarity and ease of use. The DC-60 Exp with X-Insight provides outstanding ease of use with better ergonomics, easier exams, and flexible management, even beyond your expectations.



21.5" full HD monitor



13.3" ultra-slim touch screen with angle adjustment



Interactive ultrasound APP --- MedSight/MedTouch



Unique dual-wing floating arm(optional)



Multi-Gesture operation



Elevating and rotating control panel



# May the X-Insight be with you!

As an all-around partner to trust, DC-60 Exp with X-Insight will provide you 'high efficiency with precision imaging' from all-aspect of your daily clinic application.