

## Sonosite ST

Procedural partner | Designed for efficiency



# The art of care.

Sonosite ST — envision your course of action with the clarity of a point-of-care ultrasound system designed to move with you from diagnosis to procedure.



# Make every procedure a masterpiece.

The assurance you're looking for in medical imaging is here. With its 21" all-touchscreen display, Sonosite ST is engineered to enhance the way you approach point-of-care ultrasound procedures.

Featuring an intuitive user interface for optimal workflows, Sonosite ST connects clinicians to the needs of patients with nimble performance and precision, helping you move forward with confidence.



From decades of industry insight and innovation, Sonosite ST is meticulously designed for the demands of the job. Engineered to advance from diagnosis to procedure enabling clinicians to act with confidence at every turn.



## Sonosite image clarity.

With its sleek 21" all-touchscreen interface, Sonosite ST incorporates the image clarity you have come to know and rely on from a Sonosite system.

## Workflow simplified.

Information at your fingertips — with its intuitive interface, Sonosite ST exists to guide and be guided all while improving workflows, giving you time to do more of what you do best, providing patient care.



## Streamlined design.

Sonosite ST is designed to accommodate how you move and where you go — engineered to adjust to the needs of the patient and procedure, putting all the tools in all the right places.

# Insight in hand.

## 21" portrait display

Sleek 21" all touchscreen display with large image area.

## Optimised exam types

Most-needed annotations and calculations for workflow efficiency.

## Auto Steep Needle Profiling (SNP)

In-plane needle visualisation from multiple angles designed to provide enhanced visualisation with a single touch.

## On-board education

More than 150 scan-along tutorials are on board including Anaesthesiology, Pain Management, and Musculoskeletal imaging.



## Power cord management

Prevent tangles with power cords confined within the stand column.

## Versatile storage solutions

Choose different sized storage bins and left or right-sided gel/wipes holder configurations to best suit your workflow needs.

## Optional large storage bins

Designed to keep your procedural tools organized.

## One-step adjustable height

System height can be raised or lowered for a wide range of users.

## Low-profile wheelbase

Enables you to position the system close to the bedside for better access to care.




## Data management and connectivity

Wi-Fi, ethernet port, 1080p digital video output compatible with 1080p displays, (2) USB 2.0 ports, (2) USB 3.0 ports, (1) dedicated USB port for optional printer.



## Guide needles precisely

Auto Steep Needle Profiling helps in-plane needle visualisation from multiple angles with one system interaction. The clinician simply advances the needle at the angle optimal to reach the target.

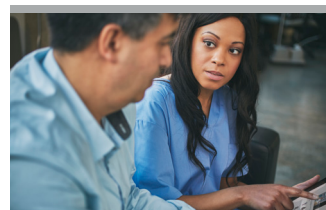

## Education at your fingertips

Sonosite is dedicated to clinician education. More than 150 scan-along education tutorials are onboard for COVID-19, diagnosis, procedures, pathology, system and transducer use.



## Rest assured with data security

Sonosite ST includes robust security features to keep patient information safe and networks secure. The system meets FIPS199, FIPS200, and NIST 800-53.

# Equipped for your procedures.



## P5-1

<b>MHz</b>	<b>Scan depth</b>
5-1 MHz Phased	5-35 cm
<b>Exam types</b>	
Abdomen	Obstetrics
Cardiac	Orbital
Focused Cardiac	Transcranial Doppler (TCD)
Lung	



## L19-5 ●●

<b>MHz</b>	<b>Scan depth</b>
19-5 MHz Linear	1-6 cm
<b>Exam types</b>	
Arterial	PIV
Lung	Superficial
Musculoskeletal	Venous
Nerve	
Ophthalmic	



## L15-4 ●

<b>MHz</b>	<b>Scan depth</b>
15-4 MHz Linear	1.5-6 cm
<b>Exam types</b>	
Arterial	Nerve
Breast	Superficial
Carotid	Venous
Musculoskeletal	



## L12-3 ●

<b>MHz</b>	<b>Scan depth</b>
12-3 MHz Linear	1.5-9 cm
<b>Exam types</b>	
Arterial	Nerve
Breast	Ophthalmic
Carotid	PIV
Lung	Superficial
Musculoskeletal	Venous



## C5-1 ●

<b>MHz</b>	<b>Scan depth</b>
5-1 MHz Curved	4.7-30 cm
<b>Exam types</b>	
Abdomen	Musculoskeletal
Gynaecology	Early Obstetrics
Lung	Obstetrics
Nerve	Spine



## C10-3 ●

<b>MHz</b>	<b>Scan depth</b>
10-3 MHz Curved	3-18 cm
<b>Exam types</b>	
Abdomen	Musculoskeletal
Gynaecology	Spine
Lung	
Nerve	



## IC10-3 ●

<b>MHz</b>	<b>Scan depth</b>
10-3 MHz Curved	3-15 cm
<b>Exam types</b>	
Gynaecology	
Early Obstetrics	
Obstetrics	
Prostate	

### Optional Needle Guide Starter Kit

- CIVCO Infiniti Plus™ In-Plane Ultrasound Needle Guide
- CIVCO AccuSITE™ Out-of-Plane Ultrasound Needle Guide
- CIVCO Disposable Endocavity Needle Guide

# Enhanced point-of-care ultrasound for your medical specialty

## Anaesthesia

### Equipping care providers with procedural insight.

Sonosite ST provides a large portrait touchscreen display, giving clinicians an immersive experience to help guide procedures and assessments such as regional nerve blocks and gastric volume.

#### **Expanded visibility.**

Sonosite image clarity helps you distinguish anatomical details, supporting needle precision and accuracy, housed in a seamless form factor, sleek panel 21" touchscreen display that's visible across the bedside.

#### **Auto Steep Needle Profiling (SNP).**

With the touch of a button, the Auto Steep Needle Profiling helps guide needle placement for more precise blocks and injections.

## Musculoskeletal

### Designed to bring clinicians closer.

Supporting Sports Medicine, Orthopaedics, Physiatry, Rheumatology and other related specialties to help you accurately diagnose and treat a variety of functional conditions involving tendons, ligaments, muscles, bursae, joints, cartilage, or peripheral nerves.

#### **Consolidated footprint.**

The sleek design and small footprint allows for easier access to patients and procedures, with effortless manoeuvring around the clinic, and features one-step height adjustment for a range of users.

#### **Immersive viewing experience.**

Sonosite image clarity on a large portrait touchscreen display allows clinicians to perform soft-tissue assessments and communicate with the patient to demystify complications and treatment needs.

## Vascular Access

### Raising the bar for procedures.

When accuracy is essential, Sonosite ST aids clinicians in placing difficult peripheral IVs and central lines, helping them guide needle placement during procedures.

#### **Workflow Simplified.**

With an intuitive user interface including integrated optimised exam types allowing clinicians to guide needle placement, easy system height adjustments, versatile storage, and a large easy-to-clean display — Sonosite ST is a proficient procedural partner.

#### **Easy to use. Easier to learn.**

A library of FUJIFILM Sonosite onboarding tools and scan-along tutorials await clinicians who utilise Sonosite ST for your ultrasound-training needs.

# Sonosite ST Ultrasound system specifications

General system and stand specifications	
Dimensions width x height unless otherwise stated	
<b>Stand footprint</b>	23 in x 25.4 in 58.4 cm x 64.5 cm
<b>Maximum height</b> (system and stand)	68 in/ 172.6 cm clinical monitor up
<b>Minimum height</b> (system and stand)	53.75 in/ 136.5 cm clinical monitor down in transport position
<b>Height adjustment range</b>	9.8 in/ 25 cm
<b>Clinical Monitor Tilt</b>	0 to 60 degrees (0 degrees tilting forward from vertical, 60 degrees tilting back from vertical)
<b>Weight</b>	109.3 lbs/ 49.6 kg
<b>System controls</b>	Touch Screen
<b>Clinical monitor</b>	
Monitor size	21.3 in/ 54.1 cm
Display size	17.01 in x 12.76 in
Image size	43.2 cm x 32.4 cm 10.03 in x 7.53 in 25.49 cm x 19.12 cm
<b>Clinical monitor resolution</b> (minimum)	
Display size	1600 x 1200 pixels
Image size	944 x 708 pixels
<b>Transducer connections</b>	(2) Transducer ports
<b>Immersion rating</b> (transducers)	IPX7
<b>Operating system</b>	Windows 10 IoT Enterprise 2019 LTSC
<b>Architecture</b>	Digital broadband
<b>Dynamic range</b>	Up to 183 dB
<b>Grey scale</b>	256 shades
<b>HIPAA compliance</b>	FIPS199, FIPS200, and NIST 800-53

## User interface

**Touch Screen on Clinical Monitor:** Projected Capacitive (PCAP) touchscreen

**Touch Screen Main Controls:** TGC (Time Gain Compensation), Gain Slider, Depth, Cine Slider, Save Clip, Save Image, Freeze, Colour, M-Mode, and 2D.

**Display formats for duplex imaging:** 1/3 and 2/3, 1/2 and 1/2, 2/3, and 1/3, and side by side duplex

## Image modes

2D Broadband Imaging, Tissue Harmonic Imaging, Pulse Inversion Harmonic Imaging, M-Mode, Colour Velocity Doppler (CVD), Colour Power Doppler (CPD), Colour Variance Mode (cardiac only)

## B-Mode Imaging Features

**Sonosite proprietary speckle reduction algorithm SonoADAPT Imaging Technology:**

Adjusts imaging parameters based on depth eliminating complicated manipulation of multiple controls.

**SonoMB:** Multi-beam technology to effectively reduce speckle artifact, enhance tissue contrast, and improve border detection.

## Auto steep needle profiling

(Available for these transducers and exams)

**C5-1:** Musculoskeletal, Nerve, Spine

**L19-5:** Arterial, Musculoskeletal, Nerve, PIV, Superficial, Venous

**L12-3:** Arterial, Breast, Carotid, Musculoskeletal, Nerve, PIV, Superficial, Venous

**L15-4:** Arterial, Breast, Carotid, MSK, Nerve, Superficial, Venous

**Needle guide starter kits:** C5-1, IC10-3, L19-5, L12-3, L15-4

## Basic measurements

### 2D Measurements:

Distance (10 Measurements), Curved Distance, Area, Circumference, Ellipse, Manual Trace, Volume, Target Depth, Angle

### M-Mode Measurements:

Distance, Time, Slope, Heart Rate

## Warranty

Sonosite offers an industry-leading, 5-year standard warranty on most ultrasound systems and Sonosite-manufactured transducers to ensure your equipment is protected. Learn more on sonosite.com

## Supported languages

**System user interface languages:** English, French, German, Italian, Spanish, Japanese, Simplified Chinese

## Connectivity and external data management

2 USB 2.0 ports (stand), 2 USB 3.0 Ports (system engine), 1 USB Port for optional printer Ethernet port, 1080p digital video output (compatible with 1080p displays) - stand only Storage capacity alert if free internal storage is less than 10%, study data Auto delete function IPv4 and IPv6 networking (DHCP, SLAAC, and static addressing), automatic time setting by synchronisation with network time server

## Wireless solutions

**2.4 GHz and/or 5 GHz wireless:** 802.11 (A, B, G, N, and AC networking)

**Key management protocols:** WPA/WPA2 personal, WPA/WPA2 enterprise (EAP-TLS, EAP-PEAP) Radius (EAP-TLS, EAP-PEAP)

## DICOM image management

Store, Modality Worklist, Modality Perform Procedure Steps (MPPS), Storage Commitment Basic Text and Comprehensive Structured Reports (SR), Supported SOP Classes: Verification, Ultrasound Image Storage, Ultrasound Multi-frame Image Storage, Secondary Capture Image Storage, Modality Worklist Information Model - FINL, Modality Performed Procedure Step, Storage Commitment Push Model, Basic Text SR, Comprehensive SR, DICOM USB export, Allow saving connectivity configuration via USB for easy replication or restoration, DICOM compliant, DICOM TLS Secure Transport Connection Profile support

## Data security

System meets Federal Information Processing Standards (FIPS) for moderate security control selections per FIPS199, FIPS200, and NIST 800-53 Secure boot

User and role based authentication LDAP support

**Security policies:** Open, WEP, Radius, WPA, WPA2 Complex password rules support

Federal System Use Notification support Guest mode support

DICOM TLS Secure Transport Connection Profile support

**Data encryption for data at rest and in transit:** FIPS 140-2 validated encryption algorithms

System wide white-listing

Hardened operating system and services System firewall

MDS2 located on Sonosite.com

## Stand includes

Input/output connections, dual gel holder, wipes holder, storage bin, lockable drawer, stand battery, cord wrap, power cord, small basket sets including dual gel holder and wipe holder, front gel holder, storage bin, power cord management, (or large basket sets including front gel holder, 4 large storage bins), stand battery, power cord

**Casters:** (4) inch swivel, 4 locking wheels

**FUJIFILM SONOSITE**  
Value from Innovation

### FUJIFILM Sonosite, Inc.

Worldwide Headquarters  
21919 30th Drive SE, Bothell, WA  
98021-3904  
Tel: +1 (425) 951 1200 or +1 (877) 657 8050  
Fax: +1 (425) 951 6800 [sonosite.com](http://sonosite.com)

### FUJIFILM Sonosite BV

European Headquarters  
Joop Geesinkweg 140, 1114 AB  
Amsterdam, Netherlands  
Tel: +31 20 462 0000

SONOSITE and the SONOSITE logo are registered and unregistered trademarks of FUJIFILM Sonosite, Inc. in various jurisdictions. FUJIFILM is a registered trademark of FUJIFILM Corporation in various jurisdictions. All other trademarks are the property of their respective owners.

Copyright © 2023 FUJIFILM Sonosite, Inc. All rights reserved. Subject to change.  
06/2023

### Sonosite Worldwide Offices

FUJIFILM Sonosite Australasia Pty Ltd: Australia . . . 1300 663 516  
FUJIFILM Sonosite Australasia Pty Ltd: New Zealand 0800 888 204  
FUJIFILM Sonosite Canada Inc. . . . . +1 888 554 5502  
FUJIFILM Sonosite GmbH-Germany . . . . . +49 69 80 88 40 30  
FUJIFILM Sonosite Iberica SL-Spain . . . . . +34 91 123 84 51  
FUJIFILM Sonosite India Pvt Ltd . . . . . +91 124 288 1100

FUJIFILM Sonosite Ltd-United Kingdom . . . . . +44 1462 341 151  
FUJIFILM Sonosite SARL-France . . . . . +33 182 880 702  
FUJIFILM Sonosite Asia Pacific Pte. Ltd. . . . . +65 6383 9933  
FUJIFILM (China) Investment Co, Ltd . . . . . +86 21 5010 6000  
FUJIFILM Medical Co, Ltd.-Japan . . . . . +81 3 0418 7190  
FUJIFILM Middle East (FZE) . . . . . +971 4 887 8722

Any patient. Anywhere. Anytime.