

ALPHA TRANCELL DELUXE

Pressure redistribution system

HUNTLEIGH



...with people in mind

The *Alpha Trancell® Deluxe* system offers highly effective alternating pressure redistribution in a low cost, simple to use package.

Using a clinically proven 10 minute cycle and a carefully calibrated pressure control on the pump, the *AlphaTrancell Deluxe* mattress replacement, overlay and seating systems are effective tools for the prevention of pressure ulcers and in the management of superficial tissue damage.

Mattress (Overlay)


- 20 fully alternating modular cells
- Durable polyurethane (PU) material
- Rapid deflate function
- Transport facility
- Securing straps

Mattress (Replacement)

- As above with additional integral foam underlay to provide full depth mattress replacement

Pump

- Manual pressure control with easy to use numbered settings
- Lightweight and compact
- Minimal noise and vibration levels
- Tough durable casing
- On/off switch
- Power on indicator
- Inflation wait indicator
- Audible and visual power fail alarm
- Audible and visual low pressure alarm
- Integrated battery back up for alarms
- Fewer recess areas where dirt can collect
- Hanging hooks
- Snap-fit tube connectors for secure connection
- Easy to use guide on pump

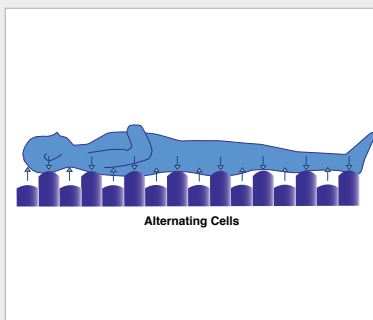

Product information				
	Mattress Overlay	Mattress Replacement	Seat Cushion	Pump
Order code	ATD002DAR	ATD001DAR	ATD003DAR	ALT01
Dimensions	2040 x 860 x 114 mm (80¼ x 34¾ x 4½")	2040 x 860 x 210 mm (80¼ x 34¾ x 8¼")	455 x 470 x 50 mm (18 x 18½ x 2")	248 x 160 x 116 mm (9¾ x 6¼ x 4½")
Weight				2.7 kg (6 lbs)
Cycle				10 minutes
Power				220-224V, 50 Hz, 14 VA
Electrical Safety				IEC60601-1–UK/Euro, IEC60601-1 – ROW
Maximum patient weight	120 kg (265 lbs)	120 kg (265 lbs)	90 kg (198 lbs)	–
Cell material	Polyurethane	Polyurethane	Polyurethane	–
Cover material	Polyurethane coated polyester	Polyurethane coated polyester	Polyurethane coated polyester	–
Warranty	2 years	2 years	1 year	2 years

Only ArjoHuntleigh designed parts, which are designed specifically for the purpose, should be used on the equipment and products supplied by ArjoHuntleigh. As our policy is one of continuous development we reserve the right to modify designs and specifications without prior notice. © and ™ are trademarks belonging to the ArjoHuntleigh group of companies. © ArjoHuntleigh, 2009

Printed in the UK • 01.TS.081.GB.11 • March 2009

HUNTLEIGH ...with people in mind

Huntleigh Healthcare UK

310-312 Dallow Road, Luton, Bedfordshire, LU1 1TD, United Kingdom

Tel: +44 (0)1582 413104

Sales: +44 (0)1582 745700

Rental: +44 (0)8457 342000 (24hr helpline)

Fax: +44 (0)1582 459100

www.huntleigh.co.uk



MEMBER OF THE GETINGE GROUP