



STERIKING[®]

IMPULSE SEALER

IS-250

Optimal choice for dental & medical clinics and laboratories

USER-FRIENDLY, ECONOMICAL

ALL-DAY READY FOR USE

CONSTANT SEALING PRESSURE, NO WARMING-UP TIME

SEALING TECHNIQUE ACCORDING TO DIN 58953, PART 7

ELECTRONIC, AUTOMATIC SEALING TIME SYSTEM

ACOUSTIC ALARM AND SIGNAL LAMP FOR THE SEALING TIME

BUILT-IN CUTTING DEVICE WITH ROTATING CIRCULAR KNIFE

PRACTICAL INTEGRATED ROLL HOLDER

SAFETY TESTED AND APPROVED BY EN, VDE

SHOCK WATCH SYSTEM FOR DELIVERY SAFETY

IMPULSE SEALER IS-250

With Integrated Roll Holder & Cutting Device


IMPULSE SEALER FOR QUICK AND RELIABLE
CLOSING OF HEAT-SEALABLE PACKAGING
MATERIALS.

SIMULTANEOUS SEALING & CUTTING

- Designed for closing all heat-sealable packaging materials in common use.
- Sealing time easily adjustable to suit the material to be sealed.
- Acoustic alarm and signal lamp indicate the adequate sealing time signal can be switched off.
- Fixed constant sealing pressure.




Sealing width
max. 250 mm

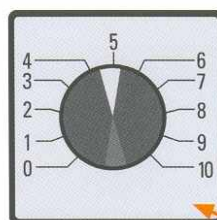

BUILT-IN CUTTING DEVICE


- Rotating circular knife enables easy manual cutting process in both directions.

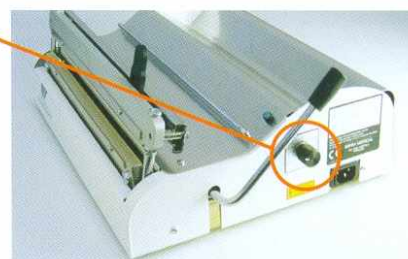

SEALING TIME SETTING

- Optimal conditions for each material can be set with the adjuster.
- Seals also Tyvek® packages well.
- Not recommended for gusseted packaging.

The IS-250 unit provides a strong seal.
Total seal width is 8 mm.

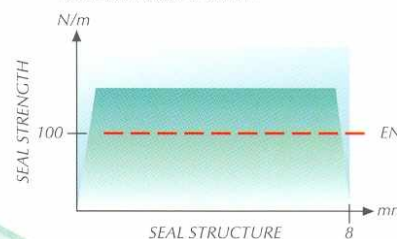

TECHNICAL DATA

Sealing time	adjustable
Seal width	8 mm, flat seal
Power consumption	870 W
Power supply	230/115 V, 50/60 Hz
Weight	10 kg
Dimensions	220 x 420 x 300 mm

WARRANTY

- 24 months from the date of purchase.

SHOCK WATCH LABEL

- Enables recognition of transport damage.
- Adhered to the transport carton.

Detailed Operating Manual and List of Spare Parts come with each unit.

WIPAK MEDICAL

Wihuri Oy Wipak
P.O. Box 45, Wipaktie 2
FI-15561 Nastola, Finland
Tel. +358 (0)20 510 311
Fax +358 (0)20 510 3333
E-Mail: info.finland@wipak.com
www.steriking.fi www.wipak.com