

Righton

Chart Projector

NP-3S



Quick and quiet

Selection of the desired chart can be completed in the blink of an eye, produces no irritating noises, and results in fast quiet examinations.

Programming

The sequence of the examinations can be programmed according to your patients' requirements, for smoother examination flow. Up to two programs can be input with 16 charts each.

Masking and red-green overlay

For either mask type, masking can be provided horizontally, vertically, or letter-by-letter. Also, all charts can be used with Red-Green overlay.



Horizontal and letter-by-letter masking also possible



Foot knobs

Four foot knobs allow flexible adjustment of projection angle in all directions

Fixed focus type

When the projector and the patient are at the same distance from the chart screen, the projector automatically produces the exact size of the chart image.

The focusing dial lets you produce sharp chart images on the screen. It's located on the bottom of the projector for easy access.



Variable focus type

Allows the placement of the projector independent of the patient. The zoom lens can be manually adjusted to obtain the necessary chart size.

Allows placement of the projector independent of the refracting distance for bright, clear images.



Standard Charts

The following four chart sets are available. These charts have been carefully selected based on importance and frequency of use, as well as space for other optional charts.

Numbers Chart 4 optional spaces

- Fixed focus type
- Semi-transparent mask
- 1.0



Europe II Chart 1 optional space

- Fixed focus type
- Solid mask
- 1.0




Snellen Numbers Chart 4 optional spaces

- Fixed and variable focus types are available.
- Semi-transparent mask
- 20/20



Multiple Chart (U.S.A.) 4 optional spaces

- Variable focus type
- Solid mask
- 20/20

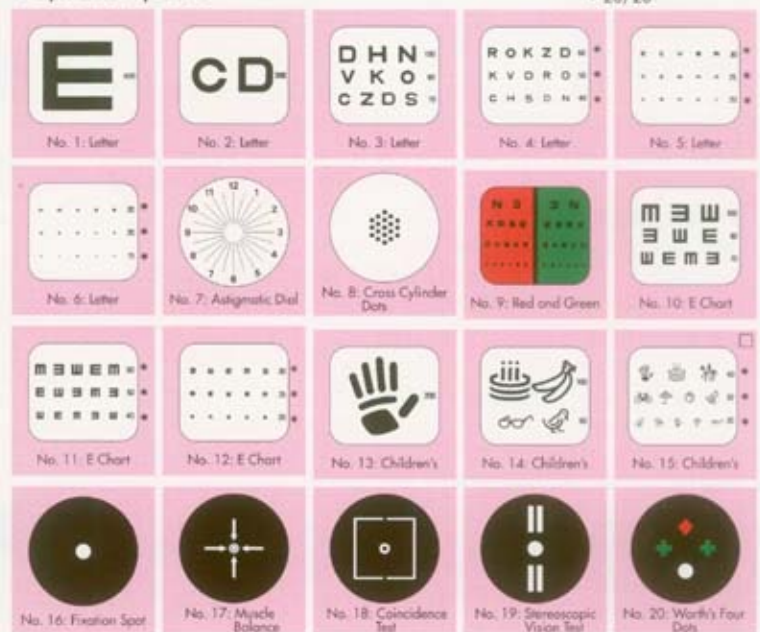

Note * Masking is possible with asterisked charts. □ Horizontal mask only.

RIGHT MFG. CO., LTD.

Ophthalmic Sales

1-47-3, Maenochi, Itabashi-ku, Tokyo 174-8633, Japan

Tel: +81-3-3960-2275 Fax: +81-3-3960-2285

e-mail: eigyousitsu@rightmfg.co.jp

