



## Less is more

**FullLife**—a minimal-contact, full-face mask—has less bulk while providing more comfort and convenience.

A revolutionary design makes FullLife a great choice for men and women who need a full-face mask but prefer something more compact with a clear line of sight.

We took a minimalist approach in developing FullLife. We made it smaller and lighter, yet ensured that it is durable. The conventional forehead pad was eliminated to provide a clear line of sight in the process.

We reduced leaks with a dual-textured, grooved cushion that creates a robust seal even with patient movement. Our straightforward, step-by-step numbered headgear system reduces application and removal effort while improving fit. With less noise, FullLife is more comfortable to live with. And with just one faceplate to fit all sizes – plus our FitPack and DuoPack options – FullLife is more comfortable for you as well. Helping you fit 100% of your patients, 100% of the time sometimes means giving you less ... and much, much more.

**PHILIPS**

**RESPIRONICS**



Elimination of the forehead pad creates a clear line of vision.

Integrated exhalation port makes operation exceptionally quiet.

One faceplate fits all 3 cushion sizes (S, M, L).

Safety entrainment valve provides quick access to room air if pressure is less than 3cm H<sub>2</sub>O.

Swivel quick-release tabs allow the hose to be easily disconnected.

Two 360° swivels allow tubing to move freely.

Branding and size indicators provide easy identification for patients and clinicians.

Chin support band ensures a secure fit and mask stability.

Cushion "ribs" add extra stability to the cushion and help prevent leaks into the eyes.


Grooved cushion creates a robust seal that accommodates movement.

SST flap and textured cushion design enhance comfort.



Step-by-step numbered headgear system makes set-up and use easy.

Dual-textured cushion surface aids assembly while preventing mask movement on the face.

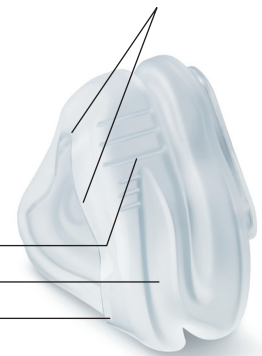

Item	Part number	U.S. HCPCS code	DuoPack (extra cushion of the same size) Part number	U.S. HCPCS code
<b>FullLife mask</b>				
Small mask with headgear	1047916	A7030,A7035	1052149	A7030,A7035,A7031
Medium mask with headgear	1047917	A7030,A7035	1052153	A7030,A7035,A7031
Large mask with headgear	1047918	A7030,A7035	1052156	A7030,A7035,A7031
Small mask without headgear	1052158	A7030	N/A	N/A
Medium mask without headgear	1052159	A7030	N/A	N/A
Large mask without headgear	1052160	A7030	N/A	N/A
FitPack: mask with headgear plus small, medium, and large cushions	1047919	A7030,A7035	N/A	N/A
<b>Replacement cushions</b>	<b>Small</b>	<b>Medium</b>	<b>Large</b>	<b>U.S. HCPCS code</b>
FullLife cushion	1047922	1047923	1047924	A7031
<b>Parts and accessories</b>			<b>Part number</b>	<b>U.S. HCPCS code</b>
FullLife headgear system			1047975	A7035

Respironics, DuoPack, FitPack, and FullLife are trademarks of Respironics, Inc. and its affiliates.



© 2009 Koninklijke Philips Electronics N.V. All rights are reserved.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

CAUTION: US federal law restricts these devices to sale by or on the order of a physician.

Hoech JJ 6/22/09 MCI 4102482 PN 1064335

Philips Healthcare is part of Royal Philips Electronics

Respironics Asia Pacific  
+65 6298 1088  
Respironics Australia  
+61 (2) 9666 4444  
Respironics Europe, Middle East, Africa  
+33 1 47 52 30 00  
Respironics United Kingdom  
+44 800 1300 845

Philips Respironics  
1010 Murry Ridge Lane  
Murrysville, PA 15668

Customer Service  
+1 724 387 4000  
800 345 6443 (toll free, US only)  
www.philips.com/respironics