

**üzümcü**






üzümcü



We are better together.  
We provide clear surgical fields for  
your accurate operations.

# SURGICAL SUCTION UNITS

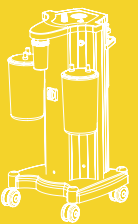

SURGICAL  
SUCTION  
UNITS



# SURGICAL SUCTION UNITS

Üzümçü Novela Series surgical suction units are life support devices immediately evacuating the blood accumulated in the surgical areas and respiratory tract by creating vacuum at a certain flow rate, increasing surgical area visibility and providing continuous respiration.



NOVELA ^



NOVELA EXTRACTOR ^

Models	Article Codes
• Novela with 2x3 L Plastic jar	33.1010
• Novela with 2x5 L Plastic jar	33.1011
• Novela with 2x3 L Glass jar	33.1012
• Novela with 2x5 L Glass jar	33.1013

Models	Article Codes
• Novela Extractor with 2x5 L Plastic Jar	33.1016

Optional Accesories	Article Codes
• Silicone Birth Cup, 50mm	33.1017
• Silicone Birth Cup, 60mm	33.1018
• Silicone Birth Cup, 70mm	33.1019
• Cannula Set, 5 pcs	33.1020

## Detailed Specifications

### Unit & Control Panel Specifications

- Rapid suction with 60 L/Min flow rate, adjustable vacuum power at 680 mmHg interval and increased surgical area visibility
- Vacuum adjustments by easily accessible control panel and vacuum indicator for monitoring vacuum level
- Maintenance free, silent, non-vibrating, oil free vacuum pump



### < Motor Protection System

- Protective technology preventing overflowing and liquid leakage to pump system by integrated buoy system and hydrophobic filter



### < Jar Specifications

- Collecting jars with scaled body, integrated easy to carry handle, 3 and 5 liter capacity options, autoclavable glass & plastic and disposable options.



### < Foot Switch Specifications


- Designed to be used in delivery operations via Malmstrom-Throe method by its integrated pedal system. Device has optional delivery assistance equipment as birth cups and fivefold cannula sets.



Fivefold Cannula Sets




Birth Cups

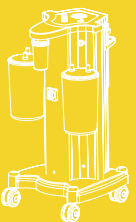

SURGICAL SUCTION UNITS

# SURGICAL SUCTION UNITS

---

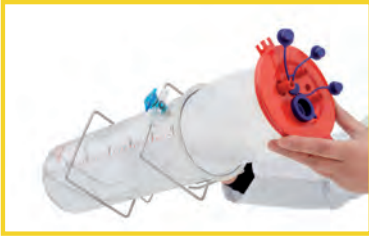

---

## < Disposable Jar

- Ease of installation and usage with optional disposable jar set



---

## < Sterile Use

- Sterilizable probe bowl to keep the suction cannula in disinfectant solution until the moment of use



---

## < Body Design

- Impact resistant and cleanable aluminium body design compatible with hospital disinfection procedures



---

## < Ergonomic Design

- Ergonomic design providing ease of carrying and moving
- Easy transfer by four castors, two of them with brake




Üzümçü HA-25 hand held suction unit is a practical device suitable for mechanical use in case of emergency operations without any electrical connection.



A.C.: 33.1024

• Max. Vacuum Power	: 450 mmHg
• Collecting Jar Capacity	: 250 ml
• Dimensions (WxLxH)	: 280 x 195 x 65mm
• Weight	: 346 gr

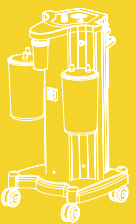

SURGICAL  
SUCTION  
UNITS

