


## Dräger ErgoStar® Consumables and Accessories

Dräger stands for modern ventilation treatment of high quality and safety when it comes to the care and treatment of patients. ErgoStar catheter mounts and airway connectors are one part in the line of breathing system equipment from Dräger.

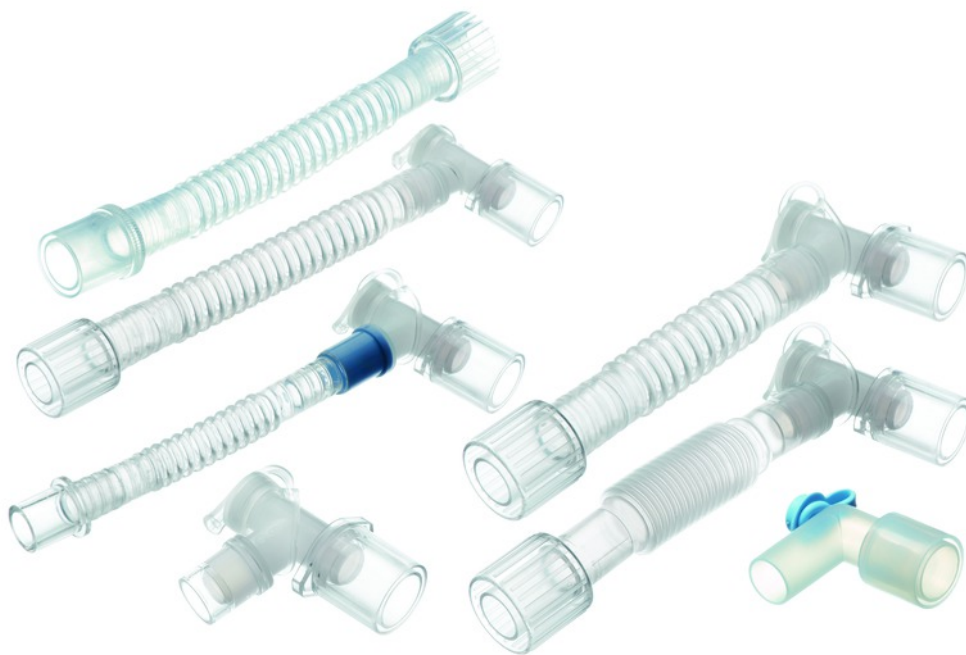

## Benefits

---

### ErgoStar

Dräger consumables are tested and manufactured to meet the highest standards in medical technology, and their functionality is convincing. Benefits of ErgoStar are:

- Smoothbore products provide low resistance to flow while external coil reduces kinks and blockages
- Extendible product provides a lightweight and flexible solution
- Double swivel connector increases ease of use
- Double cap to maintain PEEP during suctioning or bronchoscopy
- Each catheter mount and airway connector is equipped with a red safety cap so that the system is completely closed and protected until directly prior to use on the patient
- 22F/15M connector on device side as a two-in-one solution
- Standardized connectors for safe connection to other components
- Luer lock port for gas sampling
- Disposable convenience enhances workflow
- PVC-free design of gas leading components eliminates risk of softeners such as DEHP

## Technical Data

### SPECIFICATIONS

|                 | <b>ErgoStar®<br/>CM 40</b> | <b>ErgoStar®<br/>CM 45</b> | <b>ErgoStar®<br/>CM 50</b> | <b>ErgoStar®<br/>CM 55</b> | <b>ErgoStar®<br/>CM 60</b> | <b>ErgoStar®<br/>AC 90</b> | <b>ErgoStar®<br/>AC 95</b> |
|-----------------|----------------------------|----------------------------|----------------------------|----------------------------|----------------------------|----------------------------|----------------------------|
| Hose length     | 10 cm                      | 15 cm                      | 10–16 cm                   | 15 cm                      | 10 cm                      | -                          | -                          |
| Internal volume | 26 mL                      | 34 mL                      | 25–44 mL                   | 31 mL                      | 16 mL                      | 7 mL                       | 8 mL                       |
| Resistance      | 2.1 mbar at<br>60 L/min    | 2.1 mbar at<br>60 L/min    | 2.4 mbar at<br>60 L/min    | 0.9 mbar at<br>60 L/min    | 1.7 mbar at<br>30 L/min    | 2.1 mbar at<br>60 L/min    | 1.1 mbar at<br>60 L/min    |
| Leakage         | <50 mL/min at<br>60 mbar   | <50 mL/min at<br>60 mbar   | <50 mL/min at<br>60 mbar   | <25 mL/min at<br>60 mbar   | <50 mL/min at<br>60 mbar   | <50 mL/min at<br>60 mbar   | <25 mL/min at<br>60 mbar   |
| Compliance      | <1 mL/mbar at<br>60mbar    | <1 mL/mbar at<br>60mbar    | <1 mL/mbar at<br>60mbar    | <1 mL/mbar at<br>60mbar    | <1 mL/mbar at<br>60mbar    | <1 mL/mbar at<br>60mbar    | <1 mL/mbar at<br>60mbar    |
| ISO connectors  | 22M/15F –<br>22F/15M       | 22M/15F –<br>22F/15M       | 22M/15F –<br>22F/15M       | 22M/15F –<br>22F/15M       | 22M/15F –<br>15M           | 22M/15F –<br>15M           | 22M/15F –<br>15M           |
| Product         | PVC-free<br>Latex-free     | PVC-free<br>Latex-free     | PVC-free<br>Latex-free     | PVC-free<br>Latex-free     | PVC-free<br>Latex-free     | PVC-free<br>Latex-free     | PVC-free<br>Latex-free     |

Clean-room manufactured according EN ISO 14644-1:1999

## Ordering Information

|                 | <b>Packing unit (pcs.)</b> | <b>Order No.</b> |
|-----------------|----------------------------|------------------|
| ErgoStar® CM 40 | 50                         | MP01840          |
| ErgoStar® CM 45 | 50                         | MP01845          |
| ErgoStar® CM 50 | 50                         | MP01850          |
| ErgoStar® CM 55 | 50                         | MP01855          |
| ErgoStar® CM 60 | 50                         | MP01860          |
| ErgoStar® AC 90 | 100                        | MP01890          |
| ErgoStar® AC 95 | 100                        | MP01895          |

## Notes

### CORPORATE HEADQUARTERS

Drägerwerk AG & Co. KGaA  
Moislinger Allee 53–55  
23558 Lübeck, Germany  
[www.draeger.com](http://www.draeger.com)

### Manufacturer:

Dräger Medical GmbH  
Moislinger Allee 53-55  
23558 Lübeck, Germany

### As of August 2015

Dräger Medical GmbH changes  
to Drägerwerk AG & Co. KGaA

Locate your Regional  
Sales Representative at:  
[www.draeger.com/contact](http://www.draeger.com/contact)



### REGION EUROPE CENTRAL AND EUROPE NORTH

Dräger Medical GmbH  
Moislinger Allee 53-55  
23558 Lübeck, Germany  
Tel +49 451 882 0  
Fax +49 451 882 2080  
[info@draeger.com](mailto:info@draeger.com)

### REGION EUROPE SOUTH

Dräger Médical S.A.S.  
Parc de Haute Technologie  
d'Antony 2  
25, rue Georges Besse  
92182 Antony Cedex, France  
Tel +33 1 46 11 56 00  
Fax +33 1 40 96 97 20  
[d1mfr-contact@draeger.com](mailto:d1mfr-contact@draeger.com)

### REGION MIDDLE EAST, AFRICA

Dräger Medical GmbH  
Branch Office  
P.O. Box 505108  
Dubai, United Arab Emirates  
Tel +971 4 4294 600  
Fax +971 4 4294 699  
[contactuae@draeger.com](mailto:contactuae@draeger.com)

### REGION CENTRAL AND SOUTH AMERICA

Dräger Panama Comercial  
S. de R.L.  
Complejo Business Park,  
V tower, 10th floor  
Panama City  
Tel +507 377 9100  
Fax +507 377 9130  
[contactcsa@draeger.com](mailto:contactcsa@draeger.com)

### REGION ASIA / PACIFIC

Draeger Medical  
South East Asia Pte Ltd.  
25 International Business Park  
#04-27/29 German Centre  
Singapore 609916, Singapore  
Tel +65 6572 4388  
Fax +65 6572 4399  
[asia.pacific@draeger.com](mailto:asia.pacific@draeger.com)

### REGION NORTH AMERICA

Draeger Medical, Inc.  
3135 Quarry Road  
Telford, PA 18969-1042, USA  
Tel +1 215 721 5400  
Toll-free+1 800 437 2437  
Fax +1 215 723 5935  
[info.usa@draeger.com](mailto:info.usa@draeger.com)

The quality management system at Dräger Medical GmbH is certified according to  
ISO 13485, ISO 9001 and Annex II.3 of Directive 93/42/EEC (Medical devices).